


# DLITE LENS CHANGING GUIDE

## 1-10x LENS GUIDE

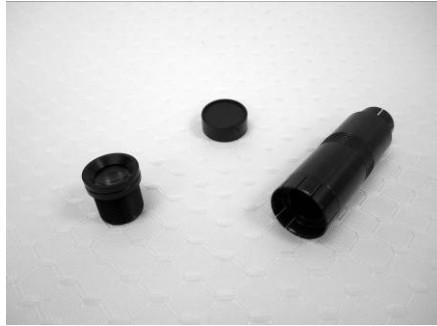

Left to Right: 1-10x Lens, Lens Cap, Lens Changing Tool

The 1-10x lens turns the Dlite into a regular video camera for use as a document camera, presentation camera, and dissecting camera. Use the instructions below to remove a lens and insert the 1-10x lens.

### Removing a lens (10-200x zoom or 500x):



1. Remove any lens tip from Dlite camera by pulling straight out.



2. Turn focus knob on the Dlite camera until lens is closest to top of camera.



3. Place small round end of blue 1-10x tool into the groove in the 10-200x lens or 4 pronged end of silver 500x tool onto 500x lens.



4. Turn the tool counter clockwise to unscrew and remove the lens.

### Inserting the 1-10x lens:



1. Remove the plastic cap from the 1-10x lens. Insert the 1-10x lens into blue 1-10x tool on large rounded side. Make sure the threaded end faces away from the tool.



2. Install the 1-10x lens by inserting and turning clockwise. **BE CAREFUL NOT TO CROSSTHREAD THE LENS.** Hold the focus knob still with thumb and tighten lens slightly with the tool. **DO NOT OVERTIGHTEN LENS.**



3. Gently pull the tool straight off the 1-10x lens. Do not twist tool while removing as this may damage the lens, pull straight off. No tip is needed when using the 1-10x lens.



4. Camera with 1-10x lens installed. Tool with lens also shown. (Use the reverse process to re-mount 10-200x lens. **DO NOT OVERTIGHTEN LENS** when re-mounting, pull tool straight off.)

1 year warranty for parts and labor. Contact us for repair or replacement within 1 year. Repair or replacement cost may apply after 1 year from invoice date.



5274 Scotts Valley Dr. #204, Scotts Valley, CA 95066  
Toll-free 877-395-1001 Fax 866-991-9591  
[www.strscopes.com](http://www.strscopes.com)